novafon®

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English

1. Key for overview	33
2. Preface	33
3. Product description	33
3.1. Scope of delivery	33
3.2. Device type, technical specifications	33
3.3. Accessories	34
3.3.1. Disc attachment	35
3.3.2. Ball attachment	35
3.3.3. Charging stand	35
3.3.4. Optional accessories	35
3.3.5. Combination with other products	35
4. Operation	36
4.1. Display	36
4.2. Initial start-up and frequency	37
4.3. Setting the intensity / strength	38
4.4. Charging	38
4.5. The App / Bluetooth	39
4.6. Attachments	40
5. Application & use & user groups	41
5.1. Intended use	41
5.2. General information for use	41
6. Contraindications and Side Effects	42
7. Technical specifications	44
7.1. Symbols used	44
7.2. Technical specifications	44
8. Important Safety Instructions	48
8.1. Warnings	48
8.2. Other precautions	52
9. Maintenance and repair	53
9.1. Cleaning and disinfection	53
9.2. Troubleshooting	54
9.2.1 Factory Settings	54
9.3. Safety check, maintenance and repair	55
9.4. Firmware updates	55
9.5. Disposal	56
10. Warranty	56
11. Customer service	57
12. Simplified EU Declaration of Conformity	57

1. Key for overview

On page 1 you will find an overview of your device with the following important elements:

- 1 Display
- 2 Multi-function button
Top: Intensity setting (choice of 3 levels)
Bottom: Frequency setting
Stand-by / ON / 100 Hz / 75 Hz / 50 Hz / OFF
- 3 Correct position of attachment on the device
- 4 Charging port

2. Preface

Congratulations on your purchase of a novafon vibration sound wave device. Our devices have delivered excellent performance and reliability for decades. Please read the following instructions carefully to ensure that you can get the most out of your novafon vibration sound wave device and enjoy using it for a long time.

Please keep these instructions for future reference.

3. Product description

3.1. Scope of delivery

Delivery includes a battery-powered novafon vibration sound wave device as well as a disc attachment, a ball attachment, a power plug, a charging cable and instructions for use.

At the Beginning of the instructions for use you will also find instructions on how to download the novafon app – for optimal use of your novafon device.

3.2. Device type, technical specifications

novafon sound wave devices have been used for decades to relieve pain and other ailments. Vibration massage can have a relaxing or invigorating effect and contribute to a sense of well-being. The transmission of gentle vibrations to the skin's surface, muscles, cartilage and tissue produces a stimulating

effect that is often still noticeable for several minutes or even hours after treatment.

Thanks to the three treatment frequencies and three intensity settings, treatment can be individually adapted to the body part and underlying symptoms in question.

The following table describes the characteristics of your new novafon device

Device type	novafon one+
Multifunction buttons	ON/ OFF (⏻) and (⏷)
frequency settings (⏻)	50 Hz, 75 Hz, 100 Hz, Standby, Off
Intensity / strength setting (⏷)	low (○) to high (⦿)
Intensity	●●●○○
Depth of penetration	6 cm +
Colour	black
Lithium ion battery	Capacity 3400 mAh, voltage 3.7 V, 12.58 Wh, Type INR18650 with protection board
Charging time	approx. 3 h
Battery life	up to 3.5 h depending on frequency and intensity

For more technical specifications please see the technical data sheet in the download section of our website.

3.3. Accessories

novafon vibration sound wave devices are equipped with two interchangeable attachments as standard.

novafon GmbH also offers other attachments that are not included in this set. Please contact your dealer or visit www.novafon.com for more information. novafon devices may only be used with original novafon accessories.

3.3.1. Disc attachment

The disc attachment is used in the treatment of larger body parts that require broader strokes. It is supplied as standard and also available as a spare part.

3.3.2. Ball attachment

The ball attachment can be used for smaller areas and targeted treatment. It is supplied as standard and also available as a spare part.

3.3.3. Charging stand

The charging stand can be purchased together with the device or separately. The charging stand offers convenient storage for your device during and after charging. These instructions for use include all the information you need on how to use and clean your stand, as well as any necessary precautionary measures.

The charging stand is delivered with an additional charging cable.

3.3.4. Optional accessories

On our website you will find additional attachments that simplify treatments in different areas of the body.

3.3.5. Combination with other products

novafon devices are not intended to be used with medical products other than novafon original accessories.

The use of accessories from other manufacturers can impair the function and performance of novafon devices.

4. Operation

4.1. Display

The display shows everything you need to know about the status and settings of your device while you are using it. Information on device status

View 1

State 1: Battery status
(only displayed at the beginning of treatment and when the battery is low)


State 2: Bluetooth status

Frequency

Runtime (current use)

Pressure warning (see view 4)

Intensity ring
(see chapter 4.2, 4.3, 4.4 , 4.5)

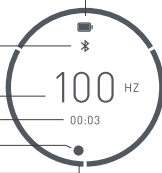View 2

Travel lock Device locked /
Switch to
unlocked (see chapter 4.1)

View 3

Welcome screen
(see chapter 4.1)
TIP: Use the app for
customisation!
Information during use

View 4

Pressure warning
(to protect the device)
(see chapter 4.3)

4.2. Initial start-up and frequency

Fully charge the device prior to switching it on for the first time (chapter 4.4). The device is delivered in locked status and must first be unlocked.

Unlock the device by pressing the ON/OFF button for 3 seconds or until the ring around the lock is completely closed. You will now see the welcome text (view 3) and then the application screen (battery level indicator, frequency, runtime; see view 1). The device starts at 100 Hz. Pressing again (○) reduces the frequency to 75 Hz, 50 Hz, or 0 Hz.

Your current frequency setting is shown on the display (see view 1).

To switch the device off, press and hold the ON / OFF button again. After about 3 seconds the device will switch off. The duration of the current treatment will be shown on the display screen before the display switches off after a further 3 seconds.

If you are carrying your novafon with you, we recommend locking the device to prevent it from accidentally switching itself on while in your bag / case. To do this, press the ON / OFF button until you see an open lock and a circle closing around the display has been completed. As soon as the circle is closed, the lock is also shown locked. The display then turns off. After 20 minutes of use, the device briefly shuts off. Devices are designed for 20 minutes of use, followed by a 15-minute break. Use longer than 20 minutes

has no medical benefit; the expected effect (pain relief, improved motor skills, etc.) is often achieved after just 5 to 10 minutes of use.

4.3. Setting the intensity / strength

You can regulate the intensity of your device by repeatedly pressing the intensity button. Here you can choose from three intensity settings: low (○), medium (○), and high (○).

You can check the intensity setting by looking at the ring around the display (see view 1). The device is set to start on the lowest intensity setting by default (○). This also corresponds to novafon's recommendations for treatment i.e. start at low intensity and increase as required, according to each individual's sensitivity.

TIP: Using the novafon app you can adapt the intensity within the individual levels.

Note: The best vibration is achieved when the device is operated without additional contact pressure. If your device automatically switches to a lower setting, or you are unable to increase the intensity, this is probably due to a low battery level. In this case, you should charge your device.

Note: If you apply a pressure force on the device during use that is potentially damaging, you will see a corresponding message on the display (see view 4). Reduce the pressure you are applying to the device to continue treatment.

4.4. Charging

To charge the device, you should only use the supplied power plug and cable, or the novafon charging stand if included (depending on the model, the charging stand may not be included in delivery). To charge your device, first connect the cable to the device or charging stand and then to the power plug

suitable for your power connection. Only then should you plug the power plug into the socket. To charge the device without a charging stand, use the charging cable included in the device packaging. Charge the device until the battery is fully charged. Do not use the novafon while it is charging. The device has a low battery indicator.

Note that only the cable supplied with the charging stand can be used to connect the charging stand. When the low battery indicator appears (see view 1), you can continue to use your novafon for a few minutes before it switches off automatically.

When using the charging stand to charge your device, make sure you are careful when placing the device onto the connector and removing it. To remove, hold the charging stand and gently release the device. The device is charging when the battery indicator appears (see view 1).

Avoid twisting the device in the charging station in order to prevent damage.

To avoid malfunctions, the device must be completely dry and free from residues of cleaning agents, disinfectants or oils / creams.

Note: When the device is fully charged, disconnect the power plug and device from the electrical supply to avoid damage. If the device is connected to the charging accessory, all pole breaking and complete disconnection from the mains can only take place if the mains plug is removed from the socket.

4.5. The App / Bluetooth

The novafon app is available in some countries. Please note that some of the medical treatment content presented there may not apply to your

device or may not be medically approved in your country.

Bluetooth is enabled by default. To respectively deactivate or activate Bluetooth, press and hold the intensity button for approx. 5 seconds. To pair your device with your mobile device and the novafon app, follow the instructions in the app in the "My NF" tab. Press the intensity and frequency buttons on the device at the same time until "Searching" appears on the display.

The color of the Bluetooth symbol (view 1) will change if the process was successful. You will now see two small triangles to the left and right of the Bluetooth symbol. To pair the device with another app or device (e.g., a second phone, tablet, etc.), the pairing process must be repeated.

The app allows you to adjust additional settings on your device:

- Customisation of display colours
- Overview of treatment durations / usage statistics
- Device information display (serial number, firmware version)
- Favourite frequency and intensity settings
- Customised welcome text on the display screen
- Firmware updates, if available, an update is only possible via the app

Note: You can also use the device without using the app. However, you will then only be able to access the basic features described in these instructions for use.

4.6. Attachments

The attachments are fixed to the device via a plug connection. They can easily be detached by turning them a quarter turn and pulling them gently off the device. The attachments click simply onto the device.

The attachments are in the correct position when the brackets on the device and the attachment lie seamlessly on top of each other. The magnets help you to find the final position. Before starting treatment, make sure that the respective attachment is engaged in the intended final position. Please refer to the drawing on page 1 for the correct position of the attachment on the device.

5. Application & use & user groups

5.1. Intended use

The novafon sound wave device is designed to create a feeling of well-being. The gentle stimulation of the vibration massage and the associated heat effect create a physiological effect that, for many users, relaxes the muscles in the contact area and can also temporarily relieve minor pain. This creates a feeling of well-being.

Note: When using the novafon one+ on sensitive parts of the body, such as the neck, face and hands, make sure to only work with low intensities.

5.2. General information for use

Treatment should not exceed a maximum of 20 minutes per session, nor a maximum of 3 sessions per day. Please note that after 20 minutes of use, the novafon should be switched off for 15 minutes in order to avoid overheating. When treating severely painful or sensitive parts of the body (e.g. face, neck, head), you should reduce / adjust the duration of treatment according to each individual's sensitivity.

Note: Application of liniments: Apply only a thin layer of creams, ointments or cosmetic / medicinal liniments to the skin, or place a fine cloth between the device and the skin. Do not allow liquids or residues to penetrate the novafon as this may cause damage. After using agents of this type, make sure you clean the device thoroughly.

Note: Discontinue treatment immediately if side effects occur (chapter 6).

The novafon appliance can be used directly on the skin or through fabric / clothing. Slowly pass the novafon over the affected area in a circular or sweeping motion. You can also use the device selectively and let it rest in place (for up to 60 seconds or until it becomes uncomfortable) to achieve a feeling of relief from pain and discomfort. Please note that you should only pass the device over the affected area with its own weight and not apply any additional pressure, as this will reduce the vibration effect. When applying to sensitive areas (e.g. face, neck, head) or multiple body parts in one session, the application time should be shortened and adjusted to individual needs. The novafon can be used on all parts of the body. Avoid contact with the temple areas or the eyes. **Do not use in the genital area.**

The session should feel comfortable and relaxing at all times. Stop using or change the position of the novafon if you find it painful or uncomfortable. **Do not use the novafon if any of the contraindications apply to you.**

It is normal for the attachments and even the device head to become warm during use. This heat is caused by the mechanics in the device. If the device head becomes hot, switch off the device and allow it to cool down for at least 15 minutes before using it again. If the device remains hot, stop using it immediately and either discard it or contact your dealer.

6. Contraindications and Side Effects

Treatments with the novafon vibration sound wave devices are not permitted for patients with:

- Open wounds /eczema (in the body region to be treated)
- Arteriosclerosis (in the body region being treated)

- Current episodes of epilepsy
- Pregnancy
- Implants (in the body region being treated)
- Heart pacemakers, defibrillators and similar
- Thromboses (in the body region being treated)
- Brain pacemakers and deep brain stimulators
- Cardiac arrhythmia
- Tumours, cancerous and not cancerous
- Acute inflammations (in the body region being treated)
- Acute episodes of inflammatory diseases

Side effects:

Please note that the following side effects may occur:

- Increased pain/unpleasant sensation (occasionally)
- Discomfort/dizziness/restlessness (occasionally)
- Skin irritation (rarely)
- Bruises (very rarely)
- Reflexes/involuntary muscle contractions (rarely)
- Muscle soreness (rarely)

Note: Pain treatment in particular can lead to an initial worsening of the symptoms before you notice any improvement.

CAUTION: Discontinue treatment immediately if you experience any side effects. Consult a physician if symptoms persist.

Note: Especially with pain management, symptoms may initially worsen before you notice improvement.

Note: Slight skin redness and warmth are not a cause for concern. This is a desired effect of local vibration therapy, caused by increased blood flow.

Do not use to treat pain or swelling of unknown origin. In this case, consult a physician for evaluation and diagnosis.

7. Technical specifications

7.1. Symbols used

- ~ Alternating current (AC)
- ≡ Direct current (DC)
- ☐ Device with protection class II according to DIN EN 61140
- Hz Hertz
- KB Short-term operation
- mA Milliampere
- mAh Milliampere-hour
- A Amp
- V Volt
- ☒ Disposal according to the EU Directive on Waste Electrical and Electronic Equipment (WEEE)
- 🏭 Displays the manufacturer of the device.
- SN Displays the serial number as defined by the manufacturer to identify a specific device.
- 🚶 Protection class of the application part: type BF
- 🌡 Temperature limit: Refers to the temperature limits to which the device can be safely exposed
- 💧 Humidity limit
- 🌬 Atmospheric pressure limit
- 📖 Refer to the instructions
- 🏠 Use inside only
- ⚠ General warning symbols
- 3 intensity settings
- ≡ Intensity button
- ⏰ Frequency setting,
Stand-by / ON / 100 Hz / 75 Hz / 50 Hz / OFF
- 📶 Bluetooth

7.2. Technical specifications

- Operating conditions: 5 °C (41 °F) to 30 °C (86 °F);
15 % to 93 % relative humidity, non-condensing;
700 hPa to 1060 hPa air pressure
- Storage / transport: -5 °C (23 °F) to 40 °C (104 °F);
15 % to 93 % relative humidity, non-condensing;
500 hPa to 1060 hPa air pressure
- Operating mode: Short-term operation 20/15 min
- Operating noise: maximum 70 dBA

Operating temperature: maximum 50 °C (122 °F) for maximum service life.

- Power plug Model 1180004-007:
- Input power: AC 100–240 V, 50–60 Hz, 300 mA
- Output power: DC 5 V, 2 A
- LOT: Shows the date of manufacture of the power supply
- A type F (EU plug) power plug is included in the standard scope of delivery. Additional mains adapters with different power plugs are optionally available or already included in the corresponding country sets.
- USB-A charging cable to spring contact pin, cable length: approx. 1 m
- One charge of a new battery guarantees at least six 20-minute treatment sessions at medium intensity. Note that frequent use at the highest intensity level can decrease battery life
- Battery: Built-in lithium-ion battery type INR18650 with protection board, life span: 800 charging cycles.
- Battery performance: Capacity 3400 mAh, voltage 3.7 V, 12.58 Wh
- IP44 protection class, which states: protected against access with objects and particles >1 mm. Protected against splashing water.

Due to the IP rating and safety requirements for devices, it is not possible for the battery to be replaced by the end user.

- Bluetooth (BLE): frequency 2.4–2.483 GHz, Modulation DSSS, modulation type GFSK 1000 kbps, effective radiated power <0.15 mW at 25 mm. Contains FCC ID: 2AC7Z-ESPMINI1 and IC: 21098-ESP32MINI1
- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and (2) this device must accept any interference

received, including interference that may cause undesired operation.

- This product complies with UK Radio Equipment Regulation 2017

This novafon device type uses RF energy for its internal function. Therefore, its RF emissions are very low and not likely to cause any interference in nearby electronic equipment.

The device type is suitable for use in all establishments, including domestic and those directly connected to the public low-voltage supply network that supplies buildings used for domestic purposes. The device may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orientating or re-locating or shielding the location. The device complies with the following immunity test levels

Immunity Test	Immunity test level
Electrostatic Discharges immunity test (ESD) IEC 61000-4-2	Contact direct discharge ±8 kV Contact indirect discharge ±8 kV Air discharge ±2 kV ±4 kV ±8 kV ±15
Power frequency magnetic field immunity test	30 A/m (50 Hz)

Immunity Test	Immunity test level
Radiated RF EM Fields and Proximity Wireless field IEC 61000-4-3	· 10/Vm from 80 MHz to 2.7 GHz 80 % AM at 1 kHz · 27 V/m, 18 Hz PM 50%, 385 MHz · 28 V/m, 18 Hz PM 50%, 450 MHz · 9 V/m, 217 Hz PM 50%, 710 MHz, 745 MHz, 780 MHz · 28 V/m, 18 Hz PM 50%, 810 MHz, 870 MHz, 930 MHz · 28 V/m, 217 Hz PM 50%, 1,720 GHz, 1,845 GHz, 1,970 GHz · 28 V/m, 217 Hz PM 50%, 2,450 GHz · 9 V/m, 217 Hz PM 50%, 5,240 GHz, 5,500 GHz, 5,785 GHz
Radiated fields in close proximity immunity test IEC 61000-4-39	· 30 kHz CW 8 A/m ; 134.2 kHz Pulse Modulation 2.1 kHz 65 A/m · 13.56 MHz Pulse Modulation 50 kHz 7.5 A/m

- Devices comply with IEC 60601-1-2:2014/AMD1:2020 and/or CISPR 11 Emission Class B and Group 1. Devices have been tested for home use in accordance with interference immunity testing levels.
- Manufacturer: NOVAFON - Elektromedizinische Geräte GmbH (novafon GmbH), Daimlerstraße 13, 71384 Weinstadt, Germany

8. Important Safety Instructions

Before starting to use the device, please read the user manual and these safety instructions carefully and strictly follow them.

8.1. Warnings

When using an electrical device basic precautions should always be followed including the following to reduce the risk of burns, fires, electric shock, possibly severe injury and possible death to the user or other people:

- ⚠ WARNING: Do not use the novafon device during charging.
- ⚠ WARNING: Always store the device in a dry and clean location.
- ⚠ WARNING: Before connecting the power plug to the mains power source, check that the mains voltage corresponds to the voltage specification on the power plug. Connect the novafon device only to a power outlet that delivers the correct voltage.
- ⚠ WARNING: To charge the device, use only the charging cable and power supply included in the delivery or the separately available charging stand. Use the correct connector plug, if applicable. Non-observance can lead to increased electromagnetic emissions or reduced electromagnetic immunity and thus to malfunctions or unexpected operational behaviour.
- ⚠ WARNING: Always unplug the device and all charging accessories from the electrical outlet after charging and before cleaning. Make sure that the mains plug is easily accessible at all times so that you can easily disconnect the device from the mains in an emergency. Do not place heavy objects or furniture in front of the power adapter, even when using the charging stand.
- ⚠ WARNING: Ensure that devices are completely dry and with no cleaning agent, disinfectant or

oil/cream residues before using the charging stand or charging cable.

- ⚠ WARNING: Do not bring the device and charging accessories (plug, charging cable, charging stand) into contact with hot surfaces (e.g. heater, hotplate, etc.). The housing of the device or adapter or the insulation of the cable could be damaged.
- ⚠ WARNING: Do not use the device when taking a shower or bath.
- ⚠ WARNING: Do not use the device where aerosol (spray) products are used, in other explosive atmospheres, or anywhere where oxygen is administered.
- ⚠ WARNING: Do not place or drop the device into water or other liquids.
- ⚠ WARNING: Do not reach for a device that has fallen into water, especially if the charging cable is plugged into a socket. Switch off the power immediately by switching off the fuse for the relevant circuit.
- ⚠ WARNING: Do not place or store the device where it could fall or be pulled into a tub, pool or sink.
- ⚠ WARNING: Protect the device, charger and charging stand from drop impact and shocks.
- ⚠ WARNING: This device may not be modified without the manufacturer's permission.
- ⚠ WARNING: novafon devices may only be used with original novafon accessories. The use of accessories from other manufacturers can impair the function and performance of novafon devices. This also applies to the power adapter and charging cables.
- ⚠ WARNING: Do not use the device and either discard or return to your dealer for examination if any of the following conditions occur:
 - if the device fails to work correctly
 - if the device is damaged, e.g. after having been dropped

- if the device was immersed in water or other liquids
- if the device exhibits other safety-related problems, such as emitting smoke, unusual or excessive heat etc.
- ⚠ WARNING: In order to prevent adverse events to the patient and operator due to electromagnetic disturbances, the novafon shall not be operated outside its intended EMC environment. Furthermore, the novafon shall not be operated if the enclosure, cables or measures for electromagnetic shielding are damaged.
- ⚠ WARNING: Do not use the power plug, charging cable or charging stand:
 - in the event of visible damage or after dropping the power plug, charging cable or charging stand
 - if the power plug, charging cable or charging stand has been exposed to extreme moisture
 - if liquid has penetrated the unit
- ⚠ WARNING: Do not charge the device for more than 24 hours at a time. Do not leave the device unattended. Disconnect the device from the power supply at night. This also applies when using the charging stand.
- ⚠ WARNING: Do not place the device on flammable surfaces when charging.
- ⚠ WARNING: Never leave the device unattended when it is turned on.
- ⚠ WARNING: Do not clamp or pinch any part of your body, especially your fingers, between the attachment and the device. During the treatment, hold the device only by the device handle.
- ⚠ WARNING: Never drop or insert any foreign object into any opening on this device. Do not insert any objects (including wires) into openings in the device.
- ⚠ WARNING: Wireless communication devices (such as Wi-Fi stations, mobile phones, cordless phones, and related base stations) can affect the device and should be kept at least 0.5 m (2 feet) away from the device during use. This also applies to cables and plug-in power supplies. Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the device, including cables specified by the Manufacturer. Otherwise, degradation of the performance of this equipment could result.
- ⚠ WARNING: Careful supervision is necessary when used on or by children, immobile or disabled individuals. novafon devices should only be used by users from the age of 13 years who have read the instructions and are competent to use the device safely.
- ⚠ WARNING: Do not use on sleeping or unconscious individuals.
- ⚠ WARNING: Keep the device out of reach of children when not in use.
- ⚠ WARNING: Small parts (e.g. cover caps) could come loose and be swallowed. Risk of suffocation!
- ⚠ WARNING: Keep children away from packaging materials. There is a risk of suffocation.
- ⚠ WARNING: Do not use on genitals, this device is not intended as an intimate pleasure toy.
- ⚠ WARNING: Do not use outdoors or in wet or damp areas.
- ⚠ WARNING: Never operate the power plug or charging stand in a high-humidity environment or in the presence of moisture. Please note that moisture is also generated by e.g. nebulisers and paraffin baths.
- ⚠ WARNING: Do not use or charge the device under blankets, pillows, towels, clothing, etc. or in places where ventilation is restricted, as it may overheat and could cause fire, electric shock or other personal injury.

- ⚠ **WARNING:** Use the device only for its intended use as described in this manual.
- ⚠ **WARNING:** Avoid contact with eyes.
- ⚠ **WARNING:** Avoid contact with mucous membranes unless using specially designed attachments.
- ⚠ **WARNING:** Do not use the device together with other devices and do not stack the devices on such products, as this may affect the performance and function of the device as well as the other product.
- ⚠ **WARNING:** Avoid overheating the device. Do not use the device for more than 20 minutes per application. Allow the device to cool down for at least 15 minutes before switching it on again.
- ⚠ **WARNING:** For technical reasons, the device can reach a temperature of up to 50 °C, especially on the head. Therefore, do not hold the device by your head for long periods of time. Change your grip position if the device feels uncomfortably warm. Switch off the device if any part of it feels hot and allow it to cool down thoroughly before using it again.

8.2. Other precautions

Note: Use the initial session to find out how your body responds to vibration therapy. If you experience any discomfort, dizziness or nausea, stop the session immediately.

Note: Do not use the device if any of the contraindications described in chapter 6 apply to you.

Note: Do not use the device for more than 20 minutes per application. Leave it to cool down for at least 15 minutes before restarting.

Note: Make sure to secure the attachments fully onto the novafon device and keep it on during treatment.

Note: If you perform multiple treatment sessions in quick succession with a maximum duration of 20 minutes and only use a single attachment, the temperature of the attachment can reach up to 50 °C. You can avoid this by removing the attachment between sessions (e.g. for cleaning) or by using alternate attachments within the session. Stop the treatment if the contact surface feels uncomfortably warm.

Note: Make sure that no substances, such as liquids, creams or the like, can penetrate into the device or the power plug.

Note: To prevent damage, make sure that the device and charging stand are standing or resting on a stable surface during the charging process.

Note: Do not twist the device in the charging stand. This may cause damage.

Note: Keep the device away from dust, lint and other particles. The penetration of dust, lint or other particles can cause the device to fail. Clean your device regularly with suitable agents to prevent any damage caused by contamination. Please note chapter 9.1.

Note: To prevent premature wear, only store the device in a dry and clean place that is protected against UV light (incl. sunlight).

Note: Use only non-alcoholic cleaning supplies. Alcohol can attack the surfaces of the housing and the attachments.

9. Maintenance and repair

9.1. Cleaning and disinfection

Do not clean the device during charging. Disconnect the device from the power supply when cleaning. If

used privately, clean the housing and attachments using a clean cloth moistened with fresh water or a mild soapy solution. For clinical and professional use, clean and disinfect device housings and all attachments with an appropriate surface disinfectant suitable for medical devices before and after each use. As a general rule, only use alcohol-free disinfectants that are suitable for use on sensitive surfaces, such as NOVACLEANER by novafon, SinAlc sensitive by HIMED or RHEOSEPT-WD plus wipes by Rheosol. If necessary, dry off with a soft cloth.

If you are using the charging stand, this can be cleaned with a clean cloth moistened with fresh water or a mild soapy solution if necessary. For professional use: To prevent contamination, only place cleaned devices in the charging stand

Ensure that devices are completely dry and with no cleaning agent, disinfectant or oil / cream residues before using the charging stand or charging cable.

9.2. Troubleshooting

If your device automatically switches to a lower setting, or you are unable to increase the intensity, this is probably due to a low battery level. In this case, you should charge your device.

You can also use the FAQs section of our website to find helpful answers to frequently asked questions. This is also available in downloadable format. If other malfunctions should occur in the device, please contact our customer service department. Never open the device yourself to rectify faults. Any attempt to repair the device yourself will void your warranty.

9.2.1 Factory Settings

You can reset your device to the factory settings by connecting it to the charging cable and then simultaneously pressing and holding the ON / OFF button

and the intensity button for approx. 10 seconds. The display will show the message "Reset". To confirm the reset, briefly press the ON / OFF button. If you do not confirm the reset, the process will be aborted.

9.3. Safety check, maintenance and repair

For the novafon device, no regular inspection or maintenance in the sense of a safety check (e.g. according to MPBetreibV) is intended. To ensure safe operation, however, you should inspect the device and attachments for obvious damage such as cracks before use. Attachments must be replaced as soon as they show clear signs of wear or damage (cracks, deformation, changes to the surface). Please note that national laws may nevertheless require you to perform safety checks on the device. Therefore, for professional use, e.g. in practices or hospitals, it is recommended to carry out an electrical safety inspection, e.g. according to DGUV V3 (German Social Accident Insurance), at least every 2 years. This generally does not apply to private individuals. If you experience problems with your device, you may dispose of it or return it to your dealer. Repairs can only be carried out by novafon or qualified dealers. Information on repair options and spare parts is available on our website. Due to cleaning requirements and the necessary IP rating, the battery cannot be replaced by the end user. The battery may only be replaced by novafon or an authorized technician.

Note: Attachments and batteries are wear parts and are therefore subject to limited warranty.

9.4. Firmware updates

Occasionally it may be necessary to update the firmware on your device. You will be informed of this either in the app or on our website. To carry out an update, your device must be paired with the app (see chapter 4.5); the app will guide you through the

necessary steps. The status of the update is displayed on the screen while it is being carried out.

9.5. Disposal

Devices must never be disposed of with common household waste. For environmentally sound disposal, the user is obligated to dispose of the old device at a designated collection point for the disposal of electrical and electronic equipment. The device contains a permanently installed battery and must be disposed of accordingly. Please observe the national legal regulations when disposing of the device. You can also send your old device to novafon; we will gladly dispose of it for you. This device is subject to EU Directive on Waste Electrical and Electronic Equipment (WEEE). To ensure the conservation of resources, please dispose of the packaging separately as waste paper.

10. Warranty

From the date of delivery, we grant a warranty of 2 years (24 months) for defects resulting from material or production defects. Warranty coverage is subject to the presentation of a proof of purchase. The battery, the attachments is guaranteed for 6 months. Excluded from the warranty: Damage due to improper use, normal wear and tear and consumption (also valid for accessories). The right of rescission (reversal of the purchase), reduction (discount on the purchase price) or subsequent delivery is initially excluded. However, you have a claim to reduction or conversion at your discretion if we do not remedy a defect that is subject to our warranty within a reasonable period of time. Any damage caused by failure to use the product in accordance with its intended use or resulting from improper or unreasonable use or maintenance is not covered by the warranty. Such damage is the sole responsibility of the user. Failure to comply with this instruction manual will void the warranty. Some states may not

allow such limitation, so it may not apply to you. This limited warranty gives you specific rights and you may have other rights which may vary from state to state. Get in touch with your dealer or supplier to claim warranty.

11. Customer service

The device may only be repaired by novafon Customer Service. Please do not send in your defective device without being asked, but rather clarify the damage and repair costs with us first. Repair could be more expensive than a new device. If the repair is no longer financially feasible within the warranty period, novafon offers a replacement unit (again with a 24-month warranty). This manual can also be viewed as a PDF file at www.novafon.com.

12. Simplified EU Declaration of Conformity

The undersigned NOVAFON – Elektromedizinische Geräte GmbH, declares that the product mentioned in this manual complies with the relevant requirements of the following EU regulations:

- Directive 2014/35/EU relating to the making available on the market of electrical equipment designed for use within certain voltage limits
- Directive 2014/30/EU relating to electromagnetic compatibility
- Directive (EU) 2015/863 amending the RoHS Directive
- Regulation (EU) 2023/1542 on batteries and waste batteries
- Regulation (EU) 2025/40 on packaging and packaging waste

The full text of the EU Declaration of Conformity is available at: <https://novafon.com/de/service/downloads#gba>

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