novafon®

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English

1. Key for overview	29
2. Preface	29
3. Product Description	29
3.1. Scope of delivery	29
3.2. Device type novafon max	29
3.3. Accessories for novafon max	30
3.3.1. Disc Attachment	31
3.3.2. Ball attachment	31
3.3.3. Charging Stand	31
3.3.4. Optional Accessories	31
3.3.5. Combination with other products	31
4. Operation	31
4.1. Display	31
4.2. Initial start-up and frequency	32
4.3. Setting the intensity/strength	33
4.4. Charging	34
4.5. The App/Bluetooth	35
4.6. Attachments	35
5. Application & use & user groups	35
5.1. Intended use	37
5.2. General informations for use	37
6. Contraindications	37
7. Technical specifications	39
7.1. Symbols used	39
7.2. Technical specifications	39
8. Important Safety Instructions	42
8.1. Warnings	42
8.2. Other precautions	47
9. Maintenance and repair	48
9.1. Cleaning and disinfection	48
9.2. Troubleshooting	49
9.2.1. Factory settings	49
9.2.2. Forced device restart (hard reset)	49
9.3. Safety check, maintenance and repair	49
9.4. Firmware updates	50
9.5. Disposal	50
10. Warranty	51
11. Customer service	51

1. Key for overview

On page 1 you will find an overview of your device

- 1 Attachment
- 2 Display
- 3 Up
- 4 On / Standby / Off
 100 Hz/ 50 Hz
- 5 Down
- 6 Charge connection

2. Preface

Congratulations on your purchase of a novafon vibration sound wave device. Our devices have delivered excellent performance and reliability for decades. Please read the following instructions carefully to ensure that you can get the most out of your novafon vibration sound wave device and enjoy using it for a long time. The explanation of the footnotes can be found on the last page of this instructions.

3. Product Description

3.1. Scope of delivery

Delivery includes a battery powered novafon max sound wave device as well as a disc attachment, a ball attachment, a power plug, a charging cable and instructions for use.

At the beginning of the instructions for use you will also find instructions on how to download the novafon App – for optimal use of your novafon device.

3.2. Device type novafon max

novafon vibration sound wave devices are designed for local vibration therapy. Different attachments are used to transmit the mechanical vibrations into the tissue. Thanks to its safe and simple application, local vibration therapy can be used in many ways.

Thanks to the treatment frequencies (50 Hz and 100 Hz) featuring five intensity settings each (low to high) treatment can be individually adapted to the body part and underlying symptoms in question.

The following table describes the characteristics of your new novafon device novafon max.

Device type	novafon max
Multifunction buttons	ON / OFF (O) and (Λ), (V)
frequency settings (O)	50 Hz, 100 Hz, Stand-by, Off
Intensity / strength setting (O)	low (O) to high (O) 5 levels
Intensity	●●●●●
Depth of penetration	9 cm+
Colour	white
Lithium ion battery	Capacity 5000 mAh, Voltage 3,6 V, 18 Wh, Type 21700 with protection board
Charging time	approx. 3.5 h
Battery life	up to 3.5 h depending on frequency and intensity

For more technical specifications please see the technical data sheet in the download section of our website.

3.3. Accessories for novafon max

novafon vibration sound wave devices are equipped with two interchangeable attachments as standard. novafon GmbH also offers other attachments that are not included in this set. Please contact your dealer or visit www.novafon.com for more information. novafon devices may only be used with original novafon accessories.

3.3.1. Disc Attachment

The disc attachment is used in the treatment of larger body parts that require broader strokes. It is supplied as standard and also available as a spare part.

3.3.2. Ball attachment

The ball attachment can be used for smaller areas and targeted treatment. It is supplied as standard and also available as a spare part.

3.3.3. Charging Stand

The charging stand can be purchased together with the device or separately. The charging stand offers convenient storage for your device during and after charging. These instructions for use include all the information you need on how to use and clean your stand, as well as any necessary precautionary measures. The charging stand is delivered with an additional charging cable.

3.3.4. Optional Accessories

On our website you will find additional attachments for various treatment goals.

3.3.5. Combination with other products

novafon devices are not intended to be used with products other than novafon original accessories.

4. Operation


4.1. Display

The display shows everything you need to know about the status

Information on device status


View 1

Battery status
Frequency
Intensity ring
(see chapter 4.2, 4.3, 4.4 , 4.5)

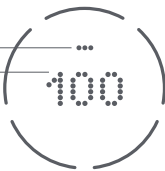View 2

Travel lock Device locked /
Switch to unlocked
(see chapter 4.2)

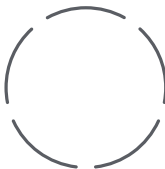View 3

Flashing rapidly: Update /
Searching for a Bluetooth device
(see Chapter 4.5)

4.2. Initial start-up and frequency

Fully charge the device prior to switching it on for the first time (chapter 4.4). The device is delivered in locked status and must first be unlocked. Unlock the device by pressing the ON/OFF button for 3 seconds or until the ring is completely open. You will now see the display (battery status, frequency, see view 1). The device starts operation at 100 Hz. Press the ON/ OFF-button again to reduce the frequency to 50 Hz or 0 Hz respectively. Your current frequency setting is shown on the display (see view 1). To switch off, press the button again for a long time, after about 3 seconds the device switches off. The display

turns off after another 3 seconds. If you are carrying your novafon with you, we recommend locking the device to prevent it from accidentally switching itself on while in your bag / case. To do this, press the ON / OFF button until the circle around the display has been completed. Then the ring lights up again. This happens with every keystroke until they completely unlock the device again.

4.3. Setting the intensity/strength

You can regulate the intensity of your device by pressing the Intensity buttons UP or DOWN. Here you can choose from five intensity settings: low (○) to high (○). You can check the intensity setting by looking at the ring around the display (see view 1). The device is set to start on the lowest intensity setting by default (○). This also corresponds to novafon's recommendations for treatment i.e. start at low intensity and increase as required, according to each individual's sensitivity.

Note: The best vibration is achieved when the device is operated without additional contact pressure. If your device automatically switches to a lower setting, or you are unable to increase the intensity, this is probably due to a low battery level. In this case, you should charge your device.

Note: In case excessive pressure is applied, the drive will lock; a click can be heard and the vibration stops. The vibration can also stop without resting by applying a lot of pressure.

If you are working exclusively at the highest intensity, the overheat protection may be activated to protect the device from damage. The device automatically switches to a lower intensity in the temperature-non-critical range. In the display, the fourth and fifth segments of the intensity ring will then flash for 5 seconds. You can then continue to use the device, but no longer at the highest intensity. As a rule, they can continue to work at low intensities, but in exceptional cases the device can be switched off completely. Note that the devices are designed for short-term

operation, i.e. after 20 minutes of use, the devices should pause for at least 15 minutes.

4.4. Charging

To charge the device, you should only use the supplied power plug and cable, or the novafon charging stand if included (depending on the model, the charging stand may not be included in delivery). To charge your device, first connect the cable to the device or charging stand and then to the power plug suitable for your power connection. Only then should you plug the power plug into the socket. Note that only the cable supplied with the charging stand can be used to connect the charging stand. To charge the device without a charging stand, use the charging cable included in the device packaging. Charge the device until the battery is fully charged (see battery icon, view 1). Do not use the novafon while it is charging. The device has a low battery indicator. When the low battery indicator appears (see view 1), you can continue to use your novafon for a few minutes before it switches off automatically. When using the charging stand to charge your device, make sure you are careful when placing the device onto the connector and removing it. To remove, hold the charging stand and gently release the device. The device is charging when the battery indicator appears (see view 1). Avoid twisting the device in the charging station in order to prevent damage. To avoid malfunctions, the device must be completely dry and free from residues of cleaning agents, disinfectants or oils/creams.

Note: When the device is fully charged, disconnect the power plug and device from the electrical supply to avoid damage. If the device is connected to the charging accessory, all pole breaking and complete disconnection from the mains can only take place if the mains plug is removed from the socket.

4.5. The App/Bluetooth

The novafon app is available in some countries. Please note that some of the medical treatment content presented there may not apply to your device or may not be medically approved in your country.

To pair your device is ready with your mobile device and the novafon app. Follow the instructions in the app using the “My NF” tab. On the device, press the intensity buttons (UP (▲) und DOWN (▼) button) at the same time until the intensity ring flashes.

Bluetooth is only activated now and only remains active as long as the device is paired and no further button is pressed. As soon as they press a button, Bluetooth is deactivated. You can use the app to perform firmware updates. The application examples shown in the app can also be carried out with the novafon max.

4.6. Attachments

The attachments are fixed to the device via a plug connection. They can easily be detached by turning them a quarter turn and pulling them gently off the device. The attachments click simply onto the device. The attachments are in the correct position when the brackets on the device and the attachment lie seamlessly on top of each other. The magnets help you to find the final position. Before starting treatment, make sure that the respective attachment is engaged in the intended final position. Due to the wave design the attachments can be individually adjusted by 90° thus allowing individual treatment direction and application to your needs. The direction of the attachment has no influence on the transmission of the vibrations into the tissue.

Please refer to the drawing on page 1 for the correct position of the attachment on the device.

5. Application & use & user groups

5.1. Intended use

The novafon sound wave device is designed to create a feeling of well-being. The gentle stimulation of the vibration massage and the associated heat effect create a

physiological effect that, for many users, relaxes the muscles in the contact area and can also temporarily relieve minor pain. This creates a feeling of well-being.

Note: When using the novafon max on sensitive parts of the body, such as the neck, face and hands, make sure to only work with low intensities.

5.2. General informations for use

Treatment should not exceed a maximum of 20 minutes per session, nor a maximum of 3 sessions per day. Please note that after 20 minutes of use, the novafon should be switched off for 15 minutes in order to avoid overheating. The device automatically switches to stand-by mode after 20 minutes of continuous use. When treating severely painful or sensitive parts of the body (e.g. face, neck, head), you should reduce / adjust the duration of treatment according to each individual's sensitivity.

Note: Application of liniments: Apply only a thin layer of creams, ointments or cosmetic / medicinal liniments to the skin, or place a fine cloth between the device and the skin. Do not allow liquids or residues to penetrate the novafon as this may cause damage. After using agents of this type, make sure you clean the device thoroughly.

Note: Discontinue treatment immediately if side effects occur (chapter 6).

The novafon appliance can be used directly on the skin or through fabric / clothing. Slowly pass the novafon over the affected area in a circular or sweeping motion. You can also use the device selectively and let it rest in place (for up to 60 seconds or until it becomes uncomfortable) to achieve a feeling of relief from pain and discomfort. Please note that you should only pass the device over the affected area with its own weight and not apply any additional pressure, as this will reduce the vibration effect. When applying to sensitive areas (e.g. face, neck, head) or multiple body parts in one session, the application time

should be shortened and adjusted to individual needs. The novafon can be used on all parts of the body. Avoid contact with the temple areas or the eyes. **Do not use in the genital area.**

The session should feel comfortable and relaxing at all times. Stop using or change the position of the novafon if you find it painful or uncomfortable. **Do not use the novafon if any of the contraindications apply to you.**

It is normal for the attachments and even the device head to become warm during use. This heat is caused by the mechanics in the device. If the device head becomes hot, switch off the device and allow it to cool down for at least 15 minutes before using it again. If the device remains hot, stop using it immediately and either discard it or contact your dealer.

6. Contraindications

Treatments with the novafon vibration sound wave devices are not permitted for patients with:

- Open wounds / eczema (in the body region to be treated)
- Arteriosclerosis (in the body region being treated)
- Current episodes of epilepsy
- Pregnancy
- Implants (in the body region being treated)
- Heart pacemakers, defibrillators and similar
- Thromboses (in the body region being treated)
- Brain pacemakers and deep brain stimulators
- Cardiac arrhythmia
- Tumours, cancerous and not cancerous
- Acute inflammations (in the body region being treated)
- Acute episodes of inflammatory diseases

Side effects:

Please note that the following side effects may occur:

- Increased pain/unpleasant sensation (occasionally)
- Discomfort/dizziness/restlessness (occasionally)
- Skin irritation (rarely)
- Bruises (very rarely)

- Reflexes/involuntary muscle contractions (rarely)
- Muscle soreness (rarely)

Note: Pain treatment in particular can lead to an initial worsening of the symptoms before you notice any improvement.

CAUTION: Discontinue treatment immediately if you experience any side effects. Consult a physician if symptoms persist.

Note: Especially with pain management, symptoms may initially worsen before you notice improvement.

Note: Slight skin redness and warmth are not a cause for concern. This is a desired effect of local vibration therapy, caused by increased blood flow.

Do not use to treat pain or swelling of unknown origin. In this case, consult a physician for evaluation and diagnosis.

7. Technical specifications

7.1. Symbols used

- ~ Alternating current (AC)
- Direct current (DC)
- ⚡ Device with protection class II according to DIN EN 61140
- Hz Hertz
- KB Short-term operation
- mA Milliampères
- mAh Milliampères per hour
- A Amp
- V Volt
- ♻ Disposal according to the EU Directive on Waste Electrical and Electronic Equipment (WEEE)
- SN Displays the serial number as defined by the manufacturer to identify a specific device.
- ⚠ Protection class of the application part: type BF

- ✂ Temperature limit: Refers to the temperature limits to which the device can be safely exposed.
- 💧 Humidity limit
- 🌬 Atmospheric pressure limit
- 📖 Refer to the instructions
- 🏠 Use inside only
- ⚠ General warning symbols
- Intensity in 5 levels
- ^ Up
- v Down
- Stand-by / ON / 100 Hz / 50 Hz / OFF

7.2. Technical specifications

- Operating conditions: 5 °C (41 °F) to 30 °C (86 °F); 15 % to 93 % relative humidity, non-condensing; 700 hPa to 1060 hPa air pressure
- Storage / transport: -5 °C (23 °F) to 40 °C (104 °F); 15 % to 93 % relative humidity, non-condensing; 500 hPa to 1060 hPa air pressure
- Operating mode: Short-term operation 20/15 min
- Operating noise: maximum 70 dBA
- Operating temperature: maximum 50 °C (122 °F) for maximum service life.
- Power plug Model 1180004-040
- Input power: AC 100–240 V, 50–60Hz, 300 mA
- Output power: DC 5 V, 2 A
- IP65 Protection class stating: Protected against dust in harmful quantities. Protection against water jets (nozzle) from any angle.
- LOT: Shows the date of manufacture of the power supply
- A type F (EU plug) power plug is included in the standard scope of delivery. Additional mains adapters with different power plugs are optionally available or already included in the corresponding country sets.
- USB-A charging cable to spring contact pin, cable length: approx. 1 m
- Battery: Built-in lithium-ion battery type 21700 with protection board, life span: 800 charging cycles. One charge of a new battery guarantees at least 6, 20-minute treat-

- ment sessions at medium intensity. Note that frequent use at the highest intensity level can decrease battery life.
- Battery performance: Capacity 5000 mAh, voltage 3,6 V, 18 Wh
- Bluetooth (BLE): frequency 2,4–2,483 GHz, Modulation DSSS, modulation type GFSK 1000kbps, effective radiated power <0.15 mW at 25 mm, Contains FCC ID: 2AC7Z-ESP-MINI1 and IC: 21098-ESP32MINI1
- This device complies with part 15 of the FCC Rules.
- Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- This product complies with UK Radio Equipment Regulation 2017
- Manufacturer: NOVAFON – Elektromedizinische Geräte GmbH (novafon GmbH), Daimlerstraße 13, 71384 Weinstadt, Germany

This novafon device type uses RF energy for its internal function. Therefore, its RF emissions are very low and not likely to cause any interference in nearby electronic equipment. The device type is suitable for use in all establishments, including domestic and those directly connected to the public low-voltage supply network that supplies buildings used for domestic purposes. The device may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orientating or re-locating or shielding the location

The device complies with the following immunity test levels

Immunity Test	Immunity test level
Electrostatic Discharges immunity test (ESD) IEC 61000-4-2	Contact direct discharge ± 8 kV Contact indirect discharge ± 8 kV Air discharge ± 2 kV ± 4 kV ± 8 kV ± 15

Immunity Test	Immunity test level
Power frequency magnetic field immunity test IEC 61000-4-8	30 A/m (50Hz)
Radiated RF EM Fields and Proximity Wire-less field IEC 61000-4-3	· 10/Vm from 80 MHz to 2.7 GHz 80 % AM at 1kHz · 27 V/m, 18 Hz PM 50%, 385 MHz · 28 V/m, 18 Hz PM 50%, 450 MHz · 9 V/m, 217 Hz PM 50%, 710 MHz, 745 MHz, 780 MHz · 28 V/m, 18 Hz PM 50%, 810 MHz, 870 MHz, 930 MHz · 28 V/m, 217 Hz PM 50%, 1,720 GHz, 1,845 GHz, 1,970 GHz · 28 V/m, 217 Hz PM 50%, 2,450 GHz · 9 V/m, 217 Hz PM 50%, 5,240 GHz, 5,500 GHz, 5,785 GHz
Radiated fields in close proximity immunity test IEC 61000-4-39	· 30 kHz CW 8 A/m ; 134.2 kHz - Pulse Modulation 2.1 kHz 65 A/m · 13.56 MHz Pulse Modulation 50 kHz 7.5 A/m

8. Important Safety Instructions

Before starting to use the device, please read the user manual and these safety instructions carefully and strictly follow them.

8.1. Warnings

When using an electrical device basic precautions should always be followed including the following to reduce the

risk of burns, fires, electric shock, possibly severe injury and possible death to the user or other people:

- ⚠ WARNING: Do not use the novafon device during charging.
- ⚠ WARNING: Always store the device in a dry and clean location.
- ⚠ WARNING: Before connecting the power plug to the mains power source, check that the mains voltage corresponds to the voltage specification on the power plug. Connect the novafon device only to a power outlet that delivers the correct voltage.
- ⚠ WARNING: To charge the device, use only the charging cable and power supply included in the delivery or the separately available charging stand. Use the correct connector plug, if applicable. Non-observance can lead to increased electromagnetic emissions or reduced electromagnetic immunity and thus to malfunctions or unexpected operational behaviour.
- ⚠ WARNING: Always unplug the device and all charging accessories from the electrical outlet after charging and before cleaning. Make sure that the mains plug is easily accessible at all times so that you can easily disconnect the device from the mains in an emergency. Do not place heavy objects or furniture in front of the power adapter, even when using the charging stand.
- ⚠ WARNING: Ensure that devices are completely dry and with no cleaning agent, disinfectant or oil/cream residues before using the charging stand or charging cable.
- ⚠ WARNING: Do not bring the device and charging accessories (plug, charging cable, charging stand) into contact with hot surfaces (e.g. heater, hotplate, etc.). The housing of the device or adapter or the insulation of the cable could be damaged.
- ⚠ WARNING: Do not use the device when taking a shower or bath.
- ⚠ WARNING: Do not use the device where aerosol (spray) products are used, in other explosive atmospheres, or anywhere where oxygen is administered.

- ⚠ WARNING: Do not place or drop the device into water or other liquids.
- ⚠ WARNING: Do not reach for a device that has fallen into water, especially if the charging cable is plugged into a socket. Switch off the power immediately by switching off the fuse for the relevant circuit.
- ⚠ WARNING: Do not place or store the device where it could fall or be pulled into a tub, pool or sink.
- ⚠ WARNING: Protect the device, charger and charging stand from drop impact and shocks.
- ⚠ WARNING: This device may not be modified without the manufacturer's permission.
- ⚠ WARNING: novafon devices may only be used with original novafon accessories. The use of accessories from other manufacturers can impair the function and performance of novafon devices. This also applies to the power adapter and charging cables.
- ⚠ WARNING: Do not use the device and either discard or return to your dealer for examination if any of the following conditions occur:
 - if the device fails to work correctly
 - if the device is damaged, e.g. after having been dropped
 - if the device was immersed in water or other liquids
 - if the device exhibits other safety-related problems, such as emitting smoke, unusual or excessive heat etc.
- ⚠ WARNING: In order to prevent adverse advents to the patient and operator due to electromagnetic disturbances, the novafon shall not be operated outside its intended EMC environment. Furthermore, the novafon shall not be operated if the enclosure, cables or measures for electromagnetic shielding are damaged.
- ⚠ WARNING: Do not use the power plug, charging cable or charging stand:
 - in the event of visible damage or after dropping the power plug, charging cable or charging stand
 - if the power plug, charging cable or charging stand has been exposed to extreme moisture
 - if liquid has penetrated the unit

- ⚠ WARNING: Never charge the device for more than 24 hours at a time. Do not leave the device unattended. Disconnect the device from the power supply at night. This also applies when using the charging stand.
- ⚠ WARNING: Do not place the device on flammable surfaces when charging.
- ⚠ WARNING: Never leave the device unattended when it is turned on.
- ⚠ WARNING: Do not clamp or pinch any part of your body, especially your fingers, between the attachment and the device. During the treatment, hold the device only by the device handle.
- ⚠ WARNING: Never drop or insert any foreign object into any opening on this device. Do not insert any objects (including wires) into openings in the device.
- ⚠ WARNING: Wireless communication devices (such as Wi-Fi stations, mobile phones, cordless phones, and related base stations) can affect the device and should be kept at least 0.5 m (2 feet) away from the device during use. This also applies to cables and plug-in power supplies. Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the device, including cables specified by the Manufacturer. Otherwise, degradation of the performance of this equipment could result.
- ⚠ WARNING: Careful supervision is necessary when used on or by children, immobile or disabled individuals. novafon devices should only be used by users from the age of 13 years who have read the instructions and are competent to use the device safely.
- ⚠ WARNING: Do not use on sleeping or unconscious individuals.
- ⚠ WARNING: Keep the device out of reach of children when not in use.
- ⚠ WARNING: Small parts (e.g. cover caps) could come loose and be swallowed. Risk of suffocation!
- ⚠ WARNING: Keep children away from packaging materials. There is a risk of suffocation.

- ⚠ WARNING: Do not use on genitals, this device is not intended as an intimate pleasure toy.
- ⚠ WARNING: Do not use outdoors or in wet or damp areas.
- ⚠ WARNING: Never operate the power plug or charging stand in a high-humidity environment or in the presence of moisture. Please note that moisture is also generated by e.g. nebulisers and paraffin baths.
- ⚠ WARNING: Do not use or charge the device under blankets, pillows, towels, clothing, etc. or in places where ventilation is restricted, as it may overheat and could cause fire, electric shock or other personal injury.
- ⚠ WARNING: Use the device only for its intended use as described in this manual.
- ⚠ WARNING: Avoid contact with eyes.
- ⚠ WARNING: Avoid contact with mucous membranes unless using specially designed attachments.
- ⚠ WARNING: Do not use the device together with other devices and do not stack the devices on such products, as this may affect the performance and function of the device as well as the other product.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.

- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. Consult the dealer or an experienced radio/TV technician for help.

- ⚠ **WARNING:** Avoid overheating the device. Do not use the device for more than 20 minutes per application. Allow the device to cool down for at least 15 minutes before switching it on again.
- ⚠ **WARNING:** For technical reasons, the device can reach a temperature of up to 50 °C, especially on the head. Therefore, do not hold the device by your head for long periods of time. Change your grip position if the device feels uncomfortably warm. Switch off the device if any part of it feels hot and allow it to cool down thoroughly before using it again.

8.2. Other precautions

Note: Use the initial session to find out how your body responds to vibration therapy. If you experience any discomfort, dizziness or nausea, stop the session immediately.

Note: Do not use the device if any of the contraindications described in chapter 6 apply to you.

Note: Do not use the device for more than 20 minutes per application. Leave it to cool down for at least 15 minutes before restarting.

Note: Make sure to secure the attachments fully onto the novafon device and keep it on during treatment.

Note: If you perform multiple treatment sessions in quick succession with a maximum duration of 20 minutes and only use a single attachment, the temperature of the attachment can reach up to 50 °C. You can avoid this by removing the attachment between sessions (e.g. for cleaning) or by using alternate attachments within the session. Stop the treatment if the contact surface feels uncomfortably warm.

Note: Make sure that no substances, such as liquids, creams or the like, can penetrate into the device or the power plug.

Note: To prevent damage, make sure that the device and charging stand are standing or resting on a stable surface during the charging process.

Note: Do not twist the device in the charging stand. This may cause damage.

Note: Keep the device away from dust, lint and other particles. The penetration of dust, lint or other particles can cause the device to fail. Clean your device regularly with suitable agents to prevent any damage caused by contamination. Please note chapter 9.1.

Note: To prevent premature wear, only store the device in a dry and clean place that is protected against UV light (incl. sunlight).

Note: Use only non-alcoholic cleaning supplies. Alcohol can attack the surfaces of the housing and the attachments.

9. Maintenance and repair

9.1. Cleaning and disinfection

Do not clean the device during charging. Disconnect the device from the power supply when cleaning. If used privately, clean the housing and attachments using a clean cloth moistened with fresh water or a mild soapy solution. If people from different households share a device, clean and disinfect device housings and all attachments with an appropriate surface disinfectant suitable for devices before and after each use. As a general rule, only use alcohol-free disinfectants that are suitable for use on sensitive surfaces, such as novacleaner by novafon, SinAlc sensitive by HIMED or RHEOSEPT-WD plus wipes by Rheosol. If necessary, dry off with a soft cloth.

If you are using the charging stand, this can be cleaned with a clean cloth moistened with fresh water or a mild soapy solution if necessary. For professional use: To prevent contamination, only place cleaned devices in the charging stand. Ensure that devices are completely dry and with no cleaning agent, disinfectant or oil/cream residues before using the charging stand or charging cable.

9.2. Troubleshooting

If your device automatically switches to a lower setting, or you are unable to increase the intensity, this is probably due to a low battery level. In this case, you should charge your device.

You can also use the FAQs section of our website to find helpful answers to frequently asked questions.

If other malfunctions should occur in the device, please contact our customer service department. Never open the device yourself to rectify faults. Any attempt to repair the device yourself will void your warranty.

9.2.1 Factory Settings

The device does not have a factory reset state.

9.2.2. Forced device restart (hard reset)

If your device no longer responds to any measures, including the factory reset described above, you can force a device restart by holding down the lower intensity button for approx. 7 seconds with the cable connected. The device will then start at 100 Hz and the lowest intensity.

9.3. Safety check, maintenance and repair

The novafon device does not require regular inspection or maintenance. To ensure safe operation, however, you should check the device regularly for obvious damage. It is therefore advisable, in the case of professional use – for example, in a doctor's surgery or hospital – to carry out an electrical safety inspection, for example in accordance with DGUV V3 (German Social Accident Insurance), at least every 2 years. Please note that there may be national laws that oblige you to carry out safety checks on the

device. This does not generally apply to private individuals. If you have any problems with your appliance either discard or return it to your dealer. Attachments shall be replaced once they show clear signs of wear or damage (cracks, deformations, change in surface). Information regarding repair options and spare parts can be found on our website.

Note: Attachments and batteries are parts subject to wear and as such are subject to a limited warranty. The battery in your device can be replaced by novafon or, where applicable, a technician authorised by novafon.

9.4. Firmware updates

Occasionally it may be necessary to update the firmware on your device. You will be informed of this either in the app or on our website. To carry out an update, your device must be paired with the app (see chapter 4.5), the app will guide you through the necessary steps. The status of the update is displayed on the screen while it is being carried out.

9.5. Disposal

Devices must never be disposed of with common household waste. For environmentally sound disposal, the user is obligated to dispose of the old device at a designated collection point for the disposal of electrical and electronic equipment. The device contains a permanently installed battery and must be disposed of accordingly. Please observe the national legal regulations when disposing of the device. You can also send your old device to novafon; we will gladly dispose of it for you. This device is subject to EU Directive on Waste Electrical and Electronic Equipment (WEEE).

To ensure the conservation of resources, please dispose of the packaging separately as waste paper.

10. Warranty

From the date of delivery, we grant a warranty of 2 years (24 months) for defects resulting from material or produc-

tion defects. Warranty coverage is subject to the presentation of a proof of purchase. The battery, charging cable, power plug and the attachments has a limited warranty of 6 months. Excluded from the warranty: Damage due to improper use, normal wear and tear and consumption (also valid for accessories). The right of rescission (reversal of the purchase), reduction (discount on the purchase price) or subsequent delivery is initially excluded. However, you have a claim to reduction or conversion at your discretion if we do not remedy a defect that is subject to our warranty within a reasonable period of time. Any damage caused by failure to use the product in accordance with its intended use or resulting from improper or unreasonable use or maintenance is not covered by the warranty. Such damage is the sole responsibility of the user. Failure to comply with this instruction manual will void the warranty. Some states may not allow such limitation, so it may not apply to you. This limited warranty gives you specific and you may have other rights which may vary from state to state.

Get in touch with your dealer or supplier to claim warranty.

11. Customer service

The device may only be repaired by novafon Customer Service. Please do not send in your defective device without being asked, but rather clarify the damage and repair costs with us first. Repair could be more expensive than a new device. If the repair is no longer financially feasible within the warranty period, novafon offers a replacement unit (again with a 24-month warranty).

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